



Soft Orthoses Directory 2020

*Allard Int. in cooperation with
Medical Specialties Inc.*



We are proud to present a small range of selected orthoses from Medical Specialties. Since 1960 they have been providing the orthopedic and athletic markets worldwide with products delivering superior performance and function. The products complement Allard Int.'s own variety of soft orthoses for a full range of options for the care of patients.



WRIST LACER™ II, Regular

A very stable wrist support that offers immobilization of the wrist with both dorsal and palmar aluminum stays that can be adjusted by medical professionals to create a custom fit. A palmar pad offers enhanced support by matching the contour of the palm.

Unique lace closure system provides equal tension across all laces to easily generate uniform compression for unparallel support. The tension is easily adjusted for patients with disproportionated wrist and forearm circumferences. The perforated suede material in combination with the polypropylene felt liner allows moisture to evaporate away from the skin and offers maximum comfort for the patient.


Length palmar side: 21 cm

Indications

Carpal Tunnel Syndrome, Wrist sprains or strains,
Post casts removal, Postoperative care.

Measurement

Wrist circumference.

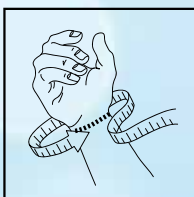

Item No. Left	Item No. Right	Color	Size	Measurement
28578 1009	28578 2009	Black	Pediatric	<12 cm
28578 1010	28578 2010	Black	X-Small	12-14 cm
28578 1011	28578 2011	Black	Small	14-16 cm
28578 1012	28578 2012	Black	Medium	16-18 cm
28578 1013	28578 2013	Black	Large	18-20 cm
28578 1014	28578 2014	Black	X-Large	20-22 cm
28578 1015	28578 2015	Black	XX-Large	22-25 cm

WRIST LACER™ II, LONG

A very stable wrist support that offers immobilization of the wrist with both dorsal and palmar aluminum stays that can be adjusted by medical professionals to create a custom fit. A palmar pad offers enhanced support by matching the contour of the palm. The long model offers even better prevention of supination/pronation of the forearm.

Unique lace closure system provides equal tension across all laces to easily generate uniform compression for unparallel support. The tension is easily adjusted for patients with disproportionated wrist and forearm circumferences. The perforated suede material in combination with the polypropylene felt liner allows moisture to evaporate and offers maximum comfort for the patient.

Length palmar side: 26 cm


Indications

Carpal Tunnel Syndrome, Wrist sprains or strains, Post casts removal, Postoperative care.

Measurement

Wrist circumference.

Item No. Left	Item No. Right	Color	Size	Measurement
28579 1011	28579 2011	Black	Small	14-16 cm
28579 1012	28579 2012	Black	Medium	16-18 cm
28579 1013	28579 2013	Black	Large	18-20 cm
28579 1014	28579 2014	Black	X-Large	20-22 cm

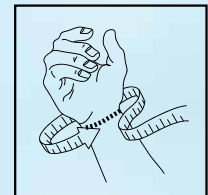

RYNO LACER™ II

A very stable wrist and thumb support with three aluminum stays - palmar, dorsal and radial, that comfortably immobilizes the wrist and thumb. Each stay is removable and can be adjusted by medical professionals to create a custom fit. Unique lace closure system provides equal tension across all laces to easily generate uniform compression for unparallel support. The tension is easily adjusted for patients with disproportionated wrist and forearm circumferences.

The perforated suede material in combination with the polypropylene felt liner allows moisture to evaporate away from the skin and offers maximum comfort for the patient. A memory foam padding in the palm and around base of the thumb enhances comfort and support.

Length palmar side: 22 cm.


Indications

DeQuervain's Tendonitis, Gamekeeper's Thumb, Postcast healing, Selected Scaphoid injuries.

Measurement

Wrist circumference.

Item No. Left	Item No. Right	Color	Size	Measurement
28576 1010	28576 2010	Black	X-Small	< 14 cm
28576 1011	28576 2011	Black	Small	14-16 cm
28576 1012	28576 2012	Black	Medium	16-18 cm
28576 1013	28576 2013	Black	Large	18-20 cm
28576 1014	28576 2014	Black	X-Large	20-22 cm

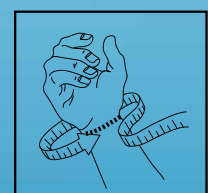



CMC-X™ Lacer Thumb Stabilizer


The Thumb Stabilizer gives a medium to high support providing pain relief for patients suffering from carpometacarpal (CMC) joint arthritis while permitting them to do the activities they enjoy. The X-Strapping System allows easily adjustment to position the first metacarpal in the desired degree of abduction to help stabilize the CMC joint while still allowing the thumb to have functional mobility. The adjustable stay pod works in concert with the “X” straps to provide gentle pressure and stabilize the CMC joint. The stay pod contains a malleable aluminum stay which is easily bent and repositioned to provide a customized fit. The lace closure system enhances ease of application. CMC-X Lacer is bilateral which reduces the carrying cost of inventory.

Indications

CMC joint osteoarthritis, Medial subluxation of the CMC joint.

Measurement

Wrist circumference.

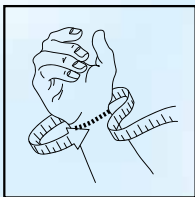

Item No.	Color	Size	L/R	Measurement
28562 0010	Black	X-Small	Bilateral	< 14 cm
28562 0011	Black	Small	Bilateral	14-16 cm
28562 0012	Black	Medium	Bilateral	16-18 cm
28562 0013	Black	Large	Bilateral	18-19,5 cm
28562 0014	Black	X-Large	Bilateral	>19,5 cm



EpiGel® Elbow Support


A comfortable, low profile design that aids in the treatment of tennis elbow (lateral epicondylitis) or golfer's elbow (medial epicondylitis). The EpiGel Elbow support provides relief with two adjustable gel pads that compress the extensor and flexor muscles without restricting circulation of the forearm. The soft gel pads also dampen vibration. A patented strap retention system allows for easy application. Unilateral design.

Indications

Lateral- and Medial Epicondylitis.

Measurement

Circumference 5 cm below the elbow.

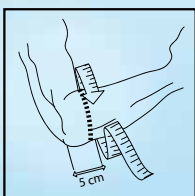

Item No.	Color	Size	Measurement
28166 0000	Black	One size	20 - 35 cm

EpiFoam™ Elbow Support

The elbow support aids in the treatment of tennis elbow (lateral epicondylitis) or golfer's elbow (medial epicondylitis). A single foam pad compresses either the extensor or flexor muscles to relieve the epicondylitis. A patented strap retention system allows for easy application. Unilateral design.


Indications

Lateral- and Medial Epicondylitis.

Measurement

Circumference 5 cm below the elbow.

Item No.	Color	Size	Measurement
28165 0000	Black	One size	20 - 35 cm

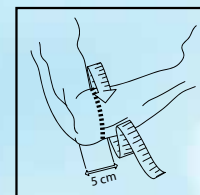

ASO® Ankle Stabilizer

The ASO® Ankle Stabilizer is a well-known, leading product in the ankle brace market worldwide. The low-profile model provides stability with stabilizing nylon straps that form a complete figure-eight to help stabilize the subtalar and talocrural joints. The straps are angled for a better contour above the malleoli for enhanced comfort and better migration resistance. A circumferential woven elastic cuff keeps laces and straps in a secure position and works in concert with the upper boot to resist the inversion or eversion forces. CoolFlex™ padding inside the boot and along the Achilles tendon offer great comfort. The shoelace is attached to the tongue to ensure lace is always centered.


Indications

Inversion or Eversion Ankle Sprains, Medial or Lateral Ankle Instability, Subtalar Joint Instability.

Measurement

Ankle circumference.

Item No.	Model	Color	Size	Measurement
28188 0009	Without stays	Black	XX-Small	22,5 - 25 cm
28188 0010	Without stays	Black	X-Small	25 - 27,5 cm
28188 0011	Without stays	Black	Small	27,5 - 30 cm
28188 0012	Without stays	Black	Medium	30 - 32,5 cm
28188 0013	Without stays	Black	Large	32,5 - 35 cm
28188 0014	Without stays	Black	X-Large	35 - 37,5 cm
28188 0015	Without stays	Black	XX-Large	37,5 - 40 cm
28188 0016	Without stays	Black	XXX-Large	> 40 cm

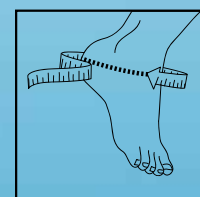



ASO® Ankle Stabilizer with Plastic Stays


The ASO® Ankle Stabilizer is a well-known, leading product in the ankle brace market worldwide. This low-profile model provides extra stability with removable plastic stays on the medial and lateral side. The stabilizing nylon straps form a complete figure-eight to help stabilize the subtalar and talocrural joints and they are angled for a better contour above the malleoli for enhanced comfort and better migration resistance. A circumferential woven elastic cuff keeps laces and straps in a secure position and works in concert with the upper boot to resist the inversion or eversion forces. CoolFlex™ padding inside the boot and along the Achilles tendon offer great comfort. The shoelace is attached to the tongue to ensure lace is always centered.

Indications

Inversion or Eversion Ankle Sprains, Medial or Lateral Ankle Instability, Subtalar Joint Instability.

Measurement

Ankle circumference.

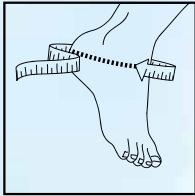

Item No.	Model	Color	Size	Measurement
28186 0010	With stays	Black	X-Small	25 - 27,5 cm
28186 0011	With stays	Black	Small	27,5 - 30 cm
28186 0012	With stays	Black	Medium	30 - 32,5 cm
28186 0013	With stays	Black	Large	32,5 - 35 cm
28186 0014	With stays	Black	X-Large	35 - 37,5 cm
28186 0015	With stays	Black	XX-Large	37,5 - 40 cm



EVO® Ankle Stabilizer


The ASO® EVO™ is a well-known, leading product in the ankle brace market worldwide. This low-profile model provides stability with stabilizing nylon straps that form a complete figure-eight to help stabilize the subtalar and talocrural joints. The straps are angled for a better contour above the malleoli for enhanced comfort and better migration resistance. The EVO model offers extra support around the calcaneus than the standard model. The circumferential woven elastic cuff is reinforced with a plastic plate to protect even better against proximal Ankle sprains. At the same time it keeps laces and straps in a secure position and works in concert with the upper boot to resist the inversion or eversion forces. CoolFlex™ padding inside the boot and along the Achilles tendon offer great comfort. The shoelace is attached to the tongue to ensure lace is always centered.

Indications

Proximal Inversion or Eversion Ankle Sprains, Medial or Lateral Ankle Instability, Subtalar Joint Instability.

Measurement

Ankle circumference.

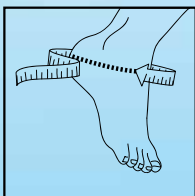

Item No.	Model	Color	Size	Measurement
28189 0010	Without stays	Black	X-Small	25 - 27,5 cm
28189 0011	Without stays	Black	Small	27,5 - 30 cm
28189 0012	Without stays	Black	Medium	30 - 32,5 cm
28189 0013	Without stays	Black	Large	32,5 - 35 cm
28189 0014	Without stays	Black	X-Large	35 - 37,5 cm
28189 0015	Without stays	Black	XX-Large	37,5 - 40 cm

Patellavator® Knee Orthosis

The Patellavator® knee orthosis has a unique design that applies pressure to the patellar tendon without creating a tourniquet around the leg or irritating the popliteal. The orthosis is able to achieve this remarkable objective by having an interlocking base strap made from Coolflex™ which is very comfortable and flexes with the knee. The pressure on the patellar tendon is applied when the external laces are pulled which force the padded rope core onto the patellar tendon. The low-profile ensures the orthosis does not irritate opposite leg. Fits easily inside pants or uniform. Unilateral.


Indications

Patellar tendonitis (Jumper's knee), Osgood Schlatter's disease, Chondromalacia of the Patella.

Measurement

Circumference below patella.

Item No.	Color	Size	Measurement
28674 0011	Black	Small	28 - 33 cm
28674 0012	Black	Medium	33 - 39 cm
28674 0013	Black	Large	39 - 44 cm
28674 0014	Black	X-Large	> 44 cm

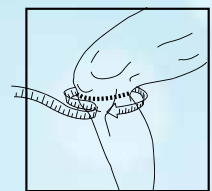

Med Spec® SI Belt

With its unique design, the SI-Belt provides exceptional pelvic-sacral support. The belt can be worn comfortably underneath clothing while remaining in its proper position throughout daily activities due to two lateral Skinloc™ Pads. Posterior plastic panels and foam pads provides additional circumferential compression to help immobilize the sacrum. Low-profile design with a Latex-free elastic liner eases the application and is comfortable enough to be worn directly against the skin.


Indications

Sacroiliac joint dysfunction, sacroiliac joint pain.

Measurement

Hip circumference.

Item No.	Color	Size	Measurement
21030 0011	Black	Small	67 - 80 cm
21030 0012	Black	Medium	80 - 98 cm
21030 0013	Black	Large	98 - 116 cm
21030 0014	Black	X-Large	116 - 133 cm
21030 0015	Black	XX-Large	133 - 150 cm

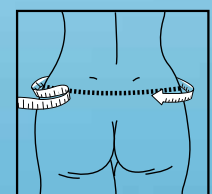

Support for Better Life!

allard | **INT**

www.allardint.com

Tel +46 42 25 27 00
Fax +46 42 25 27 25

ALLARD INTERNATIONAL
Camp Scandinavia AB
Karbingatan 38
SE-254 67 Helsingborg SWEDEN
mail.sweden@camp.se